

Vectorthotic Full Length Insoles



Features and Benefits

- This modular system brings all the advantages of custom foot orthoses without the need for casting or waiting
- Each pack contains a polypropylene composite shell, top covers and a range of medial clip-in rearfoot (2°, 4°, & 6°) and medial and lateral pre-glued forefoot (2° & 4°) posts
- The shell can be customised to the patient's requirements by clipping in the appropriate post to invert the rear foot, support the forefoot with a varus, valgus or no post as deemed suitable
- If desired, apply the moulded top cover incorporating a 2mm neoprene comfort layer, which can then be cut to the required length
- Other additions such as heel lifts or extension postings can easily be added
- Being an injection moulded span device, its strength lies in the integrity of the design and durability of the shell

Indications for Use

- Arch ache or pain
- Exercise induced kneecap pains
- Heel and ankle symptoms
- Hip and back pain on exercise
- Musculoskeletal pains linked to poor foot position and function



Ordering Information

Shoe Size (UK)	Product Code	
	Standard Top Cover (Lower arch profile)	Extra Top Cover (Higher arch profile and metatarsal support)
3 - 4	VEC05	VEC05X
5 - 6	VEC04	VEC04X
6 - 7	VEC03W	VEC03WX
7 - 8	VEC03	VEC03X
8 - 9	VEC02W	VEC02WX
9 - 10	VEC02	VEC02X
11 - 12	VEC01	VEC01X

X-Line Standard Insoles



Features and Benefits

- Features heel cradle and midfoot contoured support for maximum comfort compared to a traditional 'arch support'
- Insole is full-length and can be quickly and easily trimmed to fit most comfortable shoes
- The brushed nylon surface allows easy passage of the foot into the shoe

Indications for Use

- Arch aches or pains from foot strain
- Pain under forefoot and toe area (metatarsalgia)
- Heel and ankle pains, including mild Plantar Fasciitis (policeman's heel) and Achilles symptoms (Achilles Tendonitis)



Ordering Information

Size	Shoe Size (UK)	Product Code
Small	3 - 5	XLINE04
Medium	6 - 8	XLINE03
Large	9 - 11	XLINE02
X-Large	12 - 14	XLINE01
XX-Large	15 - 16	XLINE01X

X-Line Reinforced ³/₄



Features and Benefits

- The deep heel cup offers reassuring and effective control of the rearfoot
- There is an extra anti-pronatory feature within the heel cup so the rearfoot is actively supported
- The midfoot saddle provides essential support and improves stability throughout gait
- Unique medial arch under-bevel
- Designed to be your go-to device for many conditions particularly where footwear accommodation is an issue

Indications for Use

- Arch and foot strain
- Ankle pain
- Plantar heel pain
- Posterior Tibial Tendon Dysfunction
- Situations where full-length devices aren't a good option



Ordering Information

Size	Shoe Size (UK)	Product Code
Small	3 - 5	XLINE3/404
Medium	6 - 8	XLINE3/403
Large	9 - 11	XLINE3/402
X-Large	12 - 14	XLINE3/401

X-Line 70 Insoles



Features and Benefits

- Features heel cradle and midfoot contoured support for maximum comfort compared to a traditional 'arch support'
- Insole is full-length and can be quickly and easily trimmed to fit most comfort shoes
- The brushed nylon surface allows easy passage of the foot into the shoe
- Made from high-density foam giving a firm, highly durable foot support

Indications for Use

- Arch aches or pains from foot strain
- Pain under forefoot and toe area (metatarsalgias)
- Heel and ankle pains, including mild Plantar Fasciitis (policeman's heel) and Achilles symptoms (Achilles Tendonitis)



Ordering Information

Size	Shoe Size (UK)	Product Code
Small	3 - 5	XLINE0470
Medium	6 - 8	XLINE0370
Large	9 - 11	XLINE0270
X-Large	12 - 14	XLINE0170
XX-Large	15 - 16	XLINE01X70

X-Line Pressure Perfect Insoles



Features and Benefits

- Features special cushioning top cover material
- Heel cradle and midfoot contour supports and protects the foot by reducing heel and forefoot pressures
- Insole is full-length and can be quickly and easily trimmed to fit most comfort shoes
- Cushioning and support protects vulnerable feet from abnormal heel and forefoot pressure
- Helps keep you moving in comfort for longer

Indications for Use

- Ideal for diabetic feet
- Pain in feet due to arthritis
- Sore heels including mild Plantar Fasciitis (policeman's heel) and forefoot burning sensations (metatarsalgias)
- Active older feet
- High impact activities such as squash, badminton, netball and basketball
- Roadrunners
- Flat feet due to weight issues



Ordering Information

Size	Shoe Size (UK)	Product Code
Small	3 - 5	PRE04
Medium	6 - 8	PRE03
Large	9 - 11	PRE02
X-Large	12 - 14	PRE01
XX-Large	15 - 16	PRE01X

X-Line RIF (ReInReinforcement) Insoles



Features and Benefits

- Features plastic wrap ReInReinforcement enhancing durability and strength
- Heel cradle and midfoot contour supports the foot throughout every step providing maximum comfort compared to a traditional 'arch support'
- Insole is full-length and can be quickly and easily trimmed to fit most comfort shoes
- Strong Cambrelle® top cover for long life and enhanced insole rigidity
- Ideal for those regularly active on their feet
- Fits in any shoe with a removable in-sock/insole, most lace up shoes, trainers and walking shoes

Indications for Use

- Arch and foot strain
- Pain under forefoot and toe area (metatarsalgias)
- Ankle pains, including mild Plantar Fasciitis (policeman's heel) and Achilles symptoms
- Calf strains and other recurrent heel pain after running



Ordering Information

Size	Shoe Size (UK)	Product Code
Small	3 - 5	XLINERE04
Medium	6 - 8	XLINERE03
Large	9 - 11	XLINERE02
X-Large	12 - 14	XLINERE01
XX-Large	15 - 16	XLINERE01X

X-Line Achilles Tendonitis Insoles



Features and Benefits


- Reinforced for improved stability
- Boasts a hindfoot lift which raises the heel reducing the strain on the Achilles during the driving phase of gait
- Features a mid-foot saddle and a 'profiled arch stabilising contour', designed to reduce the speed of distal calcaneal plantarflexion of the foot during weight acceptance and prevent excess loss of navicular height

Indications for Use

- Resolves over strain and pain of the Achilles (calcaneal) tendon



X-Line TPD Insoles Firm Density

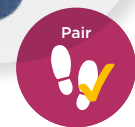

Features and Benefits

- For Tibialis Posterior Dysfunction
- Features high heel cups for heel stabilisation
- High arch support placed specifically to offload the tibialis posterior tendon
- Midfoot cradle and metatarsal support
- High plastic wrap reinforcement to strengthen insole
- Camberelle® top cover for strength and durability

Indications for Use

- Severe arch ache or pain
- Severe medial ankle and shin pain
- Any symptoms caused by large arch drop
- Perfect for immediate relief of Tibialis Posterior Dysfunction



Ordering Information

Size	Shoe Size (UK)	Product Code
Small	3 - 5	XLINEATO4
Medium	6 - 8	XLINEATO3
Large	9 - 11	XLINEATO2
X-Large	12 - 14	XLINEATO1
XX-Large	15 - 16	XLINEATO1X

Ordering Information

Size	Shoe Size (UK)	Product Code	
		Medium Density	Firm Density
Small	3 - 5	XLINETPD50NW04	XLINETPDO4
Medium	6 - 8	XLINETPD50NW03	XLINETPDO3
Large	9 - 11	XLINETPD50NW02	XLINETPDO2
X-Large	12 - 14	XLINETPD50NW01	XLINETPDO1
XX-Large	15 - 16	XLINETPD50NW01X	XLINETPDO1X

X-Line DJD Insoles



Features and Benefits

- Hallux Rigidus insoles and Hallux Limitis Insoles designed to alleviate big toe joint pain or ache during and after activity
- The first insole designed for arthritis big toe. Widely used in the NHS for the treatment of Hallux Rigidus and Hallux Limitus
- Features supports and contours that reduce the amount of bend in the big toe when walking
- Forefoot rocker precisely located to reduce big toe stress
- Heel and midfoot cradled to guarantee comfort
- The Cambrelle® top cover gives durability and an easy surface to clean

Indications for Use

- Osteoarthritis or Degenerative Joint Disease (DJD) of the big toe joint
- Big toe pain during movement
- Deep pain from bunions
- Big toe sesamoiditis
- Golf or hiking induced big toe joint pain in feet
- Excellent for use after big toe joint surgery including fusion



Ordering Information

Size	Shoe Size (UK)	Product Code
Small	3 - 5	XLINEDJD04
Medium	6 - 8	XLINEDJD03
Large	9 - 11	XLINEDJD02
X-Large	12 - 14	XLINEDJD01
XX-Large	15 - 16	XLINEDJD01X

X-Line Plantar Fasciitis Insoles



Features and Benefits

- Features a heel lift which reduces the strain on the Plantar Fascia
- Arch support to take the load from the heel
- Heel indentation and shock dot positioned to take pressure off the most sensitive area of the heel (medial process of plantar calcaneal tubercle)
- The insole also comes with an intrinsic lateral skive to the rear foot and lateral forefoot which reduces 'sudden' eversion moments
- The 1st metatarsal recess and lesser metatarsal transverse plane support the insole to improve 1st metatarsophalangeal joint function and create a stable forefoot during heel lift

Indications for Use

- To relieve heel pain caused by Plantar Fasciitis
- Relief of heel spur pain



Ordering Information

Size	Shoe Size (UK)	Product Code
Small	3 - 5	XLINEPF04
Medium	6 - 8	XLINEPF03
Large	9 - 11	XLINEPF02
X-Large	12 - 14	XLINEPF01
XX-Large	15 - 16	XLINEPF01X