

REVOLUTION SERIES™

THE INDUSTRIES MOST INNOVATIVE PROSTHETICS FOR ACTIVE AMPUTEES

Freedom Innovation is an advanced-technology prosthetics company, creating a revolution in motion as a global innovator and leader in evolving prosthetics. Merging extensive experience in composite design and prosthetic development with the industries most innovative technologies, Freedom Innovations creates the industries highest quality, top performing prosthetic products for amputees worldwide.

Runway™ Selection Guide				
Weight (lbs.)	Impact Level			Weight (kg.)
	Low	Moderate	High	
100-115	1	1	2	44-52
116-130	1	2	3	53-59
131-150	2	3	4	60-68
151-170	3	4	5	69-77
171-195	4	5	6	78-88
196-220	5	6		89-100
221-255	6			101-116

REVOLUTION SERIES™



the
Runway

The adjustable heel-height
foot fusing fashion and function



7 Studebaker, Irvine, California 92618
phone: (949) 672-0032 fax: (949) 672-0084
Order factory direct: (888) 818-6777 or through SPS (800) 767-7776
www.freedom-innovations.com

A REVOLUTION IN MOTION

A REVOLUTION IN MOTION

EXERCISE YOUR FREEDOM OF STYLE...IN ANY SHOE

Achieve the height of fashion...at any height

Runway

The Runway is the most innovative, user-adjustable heel-height foot in the industry today.

Using Freedom's highly advanced, patent-pending technology, **Runway's** simple, push-button adjustment elevates the heel up to two inches in height. The ideal combination of style and function, the **Runway** features a sandal toe, six levels of stiffness and the industries most sophisticated composite design.



Active to elegant lifestyle in one foot

The Runway

- Offers six levels of stiffness
- Incredibly smooth gait at all heel heights
- Based on the FS 2000 – one of the industries most advanced feet
- Push button adjustable heel – no tools necessary
- No hydraulic leaks or bumpers to change
- Industries highest weight limit for an adjustable heel – up to 225 lbs
- L-5981 and L-5990 reimbursements

Runway and Revolution Series are trademarks of Freedom Innovations, Inc. Model: Amy Purdy

COMPETITIVE ADVANTAGES:

The Runway's patent-pending technology:

- Offers six levels of stiffness
- Features sandal toe for versatility
- Accommodates heel heights up to 2"
- Provides full range of energy return
- Is only ankle and foot designed for the K3 (more active) lifestyle
- Features an innovative ankle that keeps the toe and heel of the foot balanced as heel height is adjusted
- Addresses cosmetic and style issues while providing maximum comfort and mobility
- Offers 18-month warranty

Freedom's Runway Foot	Competitors' Adjustable Heel Height Feet
Six levels of stiffness to provide optimal performance at all levels of weight and activity	Only one stiffness for all amputees
Designed for K-3 amputees	Accommodates K-2 amputees
Based on the design of the FS-200 foot, one of the industries most advanced prosthetic feet	Ankle attached to foot featuring dated technology
High level of energy storage and return	Minimal energy storage and return
Can accommodate weights up to 255lbs	Restricted weight limit



No heel/flat soled shoes
Lowest setting


1" heel
Medium Setting

2" heel
Highest Setting

VERSATILITY



The Runway is also compatible with your favourite pair of high heel boots.

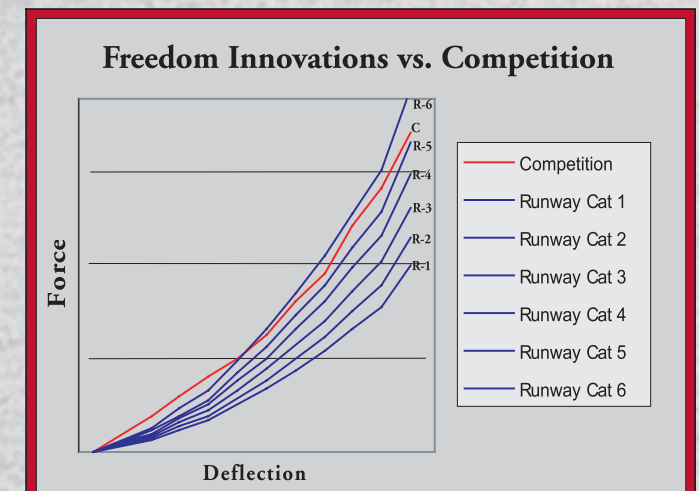

Above: Comfortable at all heel heights, the Runway does not lose performance as do competing products when worn with a high heel