



skinergyplus

 **steeper**  
Creating life's turning points, together

# Skinergy Plus Leg Covers

Manufactured using the latest innovation in silicone polymer technology, Skinergy Plus provides an extremely high level of cosmesis for the below knee user, incorporating a uniquely seamless finish.

Easy to fit, enhanced toughness and improved realism makes Skinergy Plus our most impressive leg cover yet.

19 TrueFinish™  
Colour Shades

All leg covers are available in 19 colour shades with the unique TrueFinish micropigmentation for enhanced realism.

E24654 E0 - E18 Standard colours, E24655 E0 - E18 TrueFinish colours.

Swatch colours are for indication only. To request swatch samples please contact [customerservices@steepergroup.com](mailto:customerservices@steepergroup.com)



## Durable

Integrated reinforced mesh and advanced silicone material for a more durable finish.

Significantly improved tear resistance and 25% elasticity accommodates virtually all sockets and protects the prosthetic components of the leg.

Stain resistant and easily cleaned using soap and water.

## Enhanced realism

Now available in 19 colour shades with unique TrueFinish micropigmentation for authentic realism.

Seamless finish enhances realism.

Available in both split toe and closed toe designs.

Semi-custom nail finish.

## Optimal fit for a range of users

Each Skinergy Plus comfortably stretches to fit two foot sizes and variations in foot types.

Two calf sizes available for optimal fit.

## Easy to use

Increased length by approximately 20% due to clinical feedback to provide an even smoother transition over the proximal socket brim.

Easy to use, simply roll on and off, allowing access to underlying hardware for maintenance and adjustment.




## Skinergy Plus Ordering Information

Split Toe Part Number	Closed Toe Part Number	Foot Size (cm)	Ankle Circ.		Mid Calf Circ.		Top of Cover Circ.		Cover Length From Floor to -	
			Min	Max	Min	Max	Min	Max	Mid Calf	Cover Top
TLC22LS/E*	TNS22LS/E*	22-23	20	25	28	35	26	32.5	32	55
TLC22RS/E*	TNS22RS/E*	22-23	20	25	28	35	26	32.5	32	55
TLC22LL/E*	TNS22LL/E*	22-23	20	25	32	40	30	37.5	32	55
TLC22RL/E*	TNS22RL/E*	22-23	20	25	32	40	30	37.5	32	55
TLC24LS/E*	TNS24LS/E*	24-25	22	27.5	31	39	29	36	34.5	57
TLC24RS/E*	TNS24RS/E*	24-25	22	27.5	31	39	29	36	34.5	57
TLC24LL/E*	TNS24LL/E*	24-25	22	27.5	35	44	33	41	34.5	57
TLC24RL/E*	TNS24RL/E*	24-25	22	27.5	35	44	33	41	34.5	57
TLC26LS/E*	TNS26LS/E*	26-27	22	28	32	40	30	37.5	37	61
TLC26RS/E*	TNS26RS/E*	26-27	22	28	32	40	30	37.5	37	61
TLC26LL/E*	TNS26LL/E*	26-27	22	28	37	46	35	44	37	61
TLC26RL/E*	TNS26RL/E*	26-27	22	28	37	46	35	44	37	61
TLC28LS/E*	TNS28LS/E*	28-30	24	30	34	42.5	32	40	39.5	66
TLC28RS/E*	TNS28RS/E*	28-30	24	30	34	42.5	32	40	39.5	66
TLC28LL/E*	TNS28LL/E*	28-30	24	30	39	49	37	46	39.5	66
TLC28RL/E*	TNS28RL/E*	28-30	24	30	39	49	37	46	39.5	66

Part Number	Description
NCT9740	Skinergy Plus TrueFinish 25mm Punchtool
E24655	Skinergy Plus TrueFinish Swatch Kits
E26625/E*	Skinergy Plus TrueFinish Patch Kits
E26626/E*	Skinergy Plus TrueFinish Spare Parts

**Note:** E\* denotes colour. For example, TNS22LS/E5

To access the foot bolt for maintenance purposes if required, a punch tool is available to create an access hole. The colour patch kit includes adhesive and colour matched silicone patches to provide a suitable seal. Spare patches are also available.

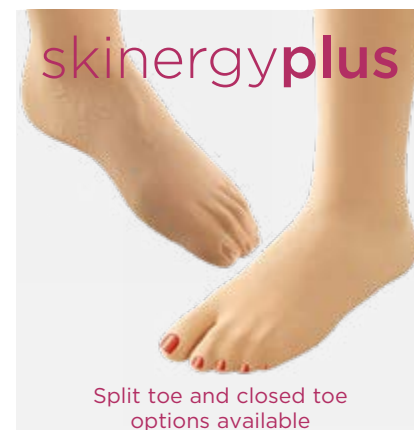

Steeper Group, Unit 3 Stourton Link, Intermezzo Drive, Leeds, LS10 1DF  
Tel: +44 (0) 870 240 4133 Email: customerservices@steepergroup.com